


■ Specifications & Features

- Maximum capacity/weight: 25kg
- Graduation: 5g
- Data comparison mark

■ Battery Installation

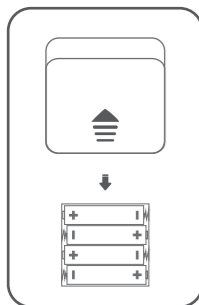

Powered by 4 x 1.5V (AAA) batteries (not included).

Remove the battery cover and dispose of packaging foil/materials before use.

Change battery: Remove the battery cover and install batteries as shown on the illustration.

Note: Please use the corresponding batteries. Remove batteries if not used for a prolonged period of time (more than 3 months).

■ Auto Off And Manual Off

For longer battery life, the scale will turn off automatically after 2 minutes without any


operation, or long press "  " to turn it off after use.

■ LCD Specification






■ Button Functions


 : Turn On/Off; ZERO

 : Short press to hold/unhold measurements, long press to switch unit.

■ Weighing

To ensure the accuracy of measurement, please use the scale on a flat and firm floor.

1. Press  button to switch on and lay down the baby carefully on the scale.
2. LCD will show baby's weight and press  to hold the measurement. The reading flashes twice for confirmation and will be stored in scale memory.
3. Comparison to previous data will be shown in difference value repeatedly. (Press  to cancel locking or allow it to turn off automatically without any operation).

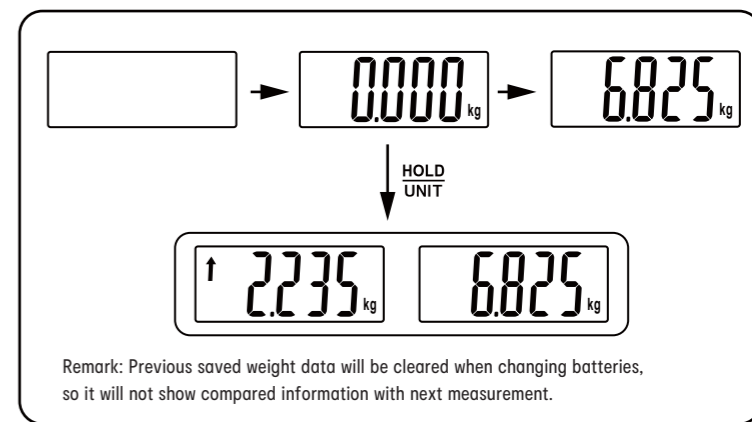

■ Indication

Warning indicator	Name	Instructions and methods
Err	OVERLOAD INDICATOR	Scale overload: the maximum capacity/weight has been reached. To avoid damage, do not overload.
Lo	LOW BATTERY INDICATOR	Battery power is low. Please replace the batteries.

■ Cleaning And Care

1. Clean the scale with a slightly damp cloth. DO NOT immerse the scale in water or use chemical/abrasive cleaning agents.
2. Always use the scale on a solid, flat surface.
3. Do not knock, impact or drop. Damage may result.
4. NOTE: This item is a precision instrument.
5. Place in a cool dry place.
6. Check the batteries if the scale cannot be turned on.
7. Designed for domestic use only.

■ CAUTION

- Do not allow baby to sleep/sit on the edge of the scale, jump on it, or bounce on the baby pan, as any result, fall or slip may cause serious injury.
- Never leave a baby unattended on the scale.
- Do not place the unit on slippery surfaces.
- Do not immerse into water or pour liquids onto the scale.

12 Month Warranty

Thank you for your purchase from Kmart.

Kmart Australia Ltd warrants your new product to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

Kmart will provide you with your choice of a refund, repair or exchange (where possible) for this product if it becomes defective within the warranty period. Kmart will bear the reasonable expense of claiming the warranty. This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre on 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or alternatively, via Customer Help at Kmart.com.au for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to our Customer Service Centre at 690 Springvale Rd, Mulgrave Vic 3170.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

KEYCODE: 42963318

MADE IN CHINA

FOR AU / NZ: IMPORTED FOR KMART STORES IN AUSTRALIA AND NEW ZEALAN

KMART AUSTRALIA - 690 SPRINGVALE ROAD, MULGRAVE, VIC 3170 AUSTRALIA.

KMART NEW ZEALAND - REGIONAL OFFICE C/O KMART PAPATOETOE STORE,
HUNTERS PLAZA, GREAT SOUTH ROAD, PAPATOETOE, AUCKLAND, NEW ZEALAND.

anko

INSTRUCTION MANUAL

DIGITAL BABY SCALE



PLEASE READ CAREFULLY BEFORE USE AND KEEP FOR FUTURE REFERENCE.

The illustrations shown in this manual are for reference only. The manufacturer reserves the right to change any specification or feature without further notice.