


Attention

- The product is not waterproof, can not be put into the water for cleaning, in order to avoid damage.
- The product can be powered by a rechargeable battery, but don't charge the battery with the product.
- If you do not use the product for a long time, please remove the battery, in order to avoid damage to battery caused by the product, then put the product in packaging box, and placed in a cool, dry and ventilated place.

Technical Specifications

Rated input voltage	DC4.5V~DC5V
Rated power	25W
Power supply	AAA battery / USB
Product specifications	81X81X90 (mm)

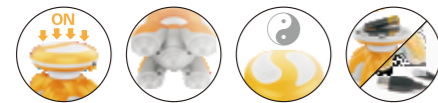

Made in China

ZhenQiang XF-68 Luminous Taiji Massager

Product Usage Instruction

Please read the instruction carefully before using and keep it well.




Product Introduction

- The product is a hand-held electric shock massager, with beautiful appearance, comfortable feel, easy to operation, suitable for massage in many parts of the body, it is a good companion for life, study and work.
- When the product is operate, the vibration effect can reach the deep level of the muscle, and achieve the massage effect of kneading, pinch, pushing and rubbing.
- If you use of the product for long term, it can free the channels and network vessels, vital energy and blood, eliminate soreness and fatigue, improve physiological functions, and enhance the resistance to disease for body
- The product uses high quality silent motor, the shell adopts environmental protection engineering material, with stable performance, health and safety.

Product Features


- The product is very easy to use, just hold the top, gently press down, it can automatically open the shock massage mode.
- It compared with the general three-leg massager, adds a massage head, to enhance more than 30% massage efficiency.
- It compared with the general mini massager, the product has better styling and comfortable hand feeling.
- The top of Taiji has massage salient points, it can massage your hand palm when in use.
- The area of Taiji salient point will show the light effect, cool luminous in the dark area.
- It support dry battery and USB line power supply.
- The product is small, easy to carry and storage.

Product parts introduction

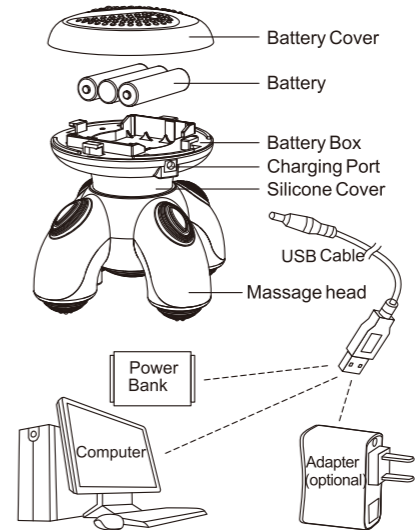


Product Usage

- Use a dry battery (AAA) for power supply

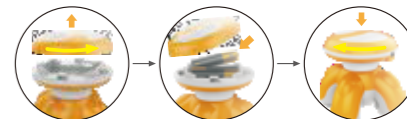

Unscrew battery cover Fit in AAA battery Screw down battery cover

- Use the USB Cable (without batteries) for power supply

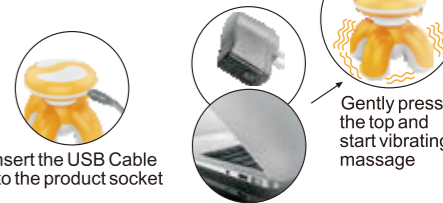

Insert the USB Cable into the product socket

Gently press the top and start vibrating massage

Insert the USB Cable into power adapter or computer or power bank