


# anko

## Cordless Shiatsu Massager

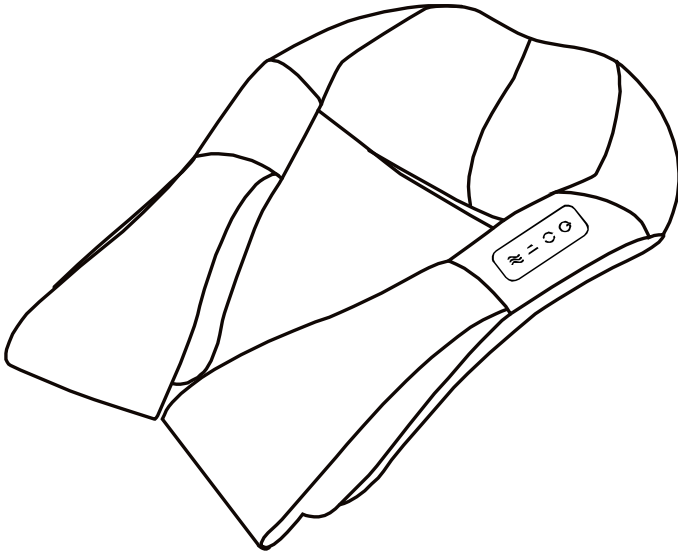

Model No.: CF-6302RNC

Rev 5

**Note:** The specifications and/or components of this appliance are subject to change without prior notice.

<b>1</b>	<b>Warnings</b> .....	<b>3</b>
<b>2</b>	<b>Important Safeguards</b> .....	<b>3</b>
<b>3</b>	<b>Introduction</b> .....	<b>6</b>
<b>4</b>	<b>Components</b> .....	<b>7</b>
<b>5</b>	<b>Operation Procedure</b> .....	<b>7</b>
<b>5.1</b>	<b>Charging</b> .....	<b>7</b>
<b>5.2</b>	<b>Massage Control Operation</b> .....	<b>8</b>
<b>6</b>	<b>How to Position your Massager</b> .....	<b>9</b>
<b>7</b>	<b>Cleaning the Massager</b> .....	<b>10</b>
<b>8</b>	<b>Completion and Storage</b> .....	<b>10</b>
<b>9</b>	<b>Technical Specifications</b> .....	<b>10</b>
<b>10</b>	<b>Troubleshooting</b> .....	<b>11</b>
<b>11</b>	<b>Warranty Against Defect</b> .....	<b>12</b>

**READ AND SAVE THESE INSTRUCTIONS FOR  
FUTURE REFERENCE**

# 1 Warnings

We recommend checking with your doctor before using the massager, especially if you have health problems.

Before using this appliance, read and follow all warnings and instructions in this manual, even if you are familiar with this product.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

# 2 Important Safeguards



Please read this instruction manual carefully before first use. It is advisable to keep this instruction manual for future reference.

**DO NOT** use this product as a substitute for medical attention.

The appliance has heated surfaces. Persons insensitive to heat must be careful when using the appliance.

Prolonged Heat application can result in skin burns.

**DO NOT** use this product for more than 15 minutes at a time. Allow at least 30 minutes cool down period in between uses. Extensive use could lead to the appliance overheating and reduce the life of the product.

The unattended use by children or incapacitated persons may be dangerous.



This product is not to be used by very young children (0-3 years)

This product is not intended for medical use. It is intended only to provide a comfortable massage.

If you have any concerns regarding your health, consult your doctor before using this appliance.

Individuals with pacemakers should consult a physician before using this product.

In case of pregnancy, diabetes or illness, consult your doctor before using this product.

Massage should be pleasant and comfortable. Should pain or discomfort result, or if you experience pain in a muscle or joint for a prolonged period of time, stop use and consult your doctor.

Never use any massager on the face, open wounds, discoloured areas, areas of poor circulation or any areas of the body that are swollen, burned, inflamed or where skin eruptions or sores are present.

**DO NOT** use the appliance while under the influence of alcohol, sleeping pills or sedatives.

It is normal to feel itching, tingling or see redness on area of the skin where the massager was used after use. These symptoms will soon disappear. However, if these symptoms remain for an extended period, we suggest you consult your doctor as there may be a more serious underlying issue.

To prevent bruising, do not press too hard on the device.

**DO NOT** leave the appliance unattended when in use or plugged in to charge.

Take care that children cannot touch the appliance when in use, and make sure they cannot reach and pull the battery charger power cord.

To protect against the risk of electric shock, do not immerse the appliance, cord or battery charger in water or other liquids.

**DO NOT** use or store the product near sinks, baths or showers containing water or near running water.

**DO NOT** use in the bath, shower, spa. Pool or other humid / damp areas etc.

Always disconnect the battery charger from the power point after use.

**DO NOT** use this appliance if it has been dropped or any part of the appliance (including cord or battery charger) has been damaged.

The use of accessories not intended for use with this appliance may cause injuries to the user or damage to the appliance.

**DO NOT** place this unit on other appliances, on uneven surfaces or where it could be subject to heat sources (e.g. radiators or stoves), direct sunlight, excessive dust

This appliance should not be placed near a hot gas or electrical burner or placed in a heated oven.

The appliance should always be switched off before disconnecting from the power supply. Never pull the cord to disconnect the battery charger from either the appliance or the power point. Grasp the plug / battery charger body.

There are no user serviceable parts in this appliance.

If this appliance is to be used by a third party, please supply the instruction manual with it.

Ensure children and babies do not play with plastic bags or any of the packaging material.

**DO NOT** operate the appliance with wet hands or wet feet to reduce the risk of electric shock.

**DO NOT** carry the appliance by the battery charger power cord.

**DO NOT** use the appliance while sleeping or if feeling sleepy.

**DO NOT** use this product before bed. The massager has a stimulating effect and may delay sleep.

**DO NOT** use this product while in bed.

**DO NOT** use machinery or heavy equipment while using this appliance.

**DO NOT** cover the appliance while in use. The product will overheat and may cause fire, electrocution, damage to property or injury to persons.

This appliance is intended for human use only. Not for animals and pets.

**DO NOT** Allow Pets or Animals to Lie or Chew on this appliance or on the battery charger electrical cord etc.

Never reach or touch an appliance which has fallen into water. Unplug from the mains power point immediately.

Position the battery charger cord properly so that it will not be walked on, tripped over or pinched by items placed on or against it.

**DO NOT** store or use the appliance near or in high temperature areas.

If the battery charger or battery charger power cord is damaged, do not use. It must be disposed in order to avoid a hazard.

**ONLY** for indoor domestic use. Industrial or commercial use invalidates the warranty.

No liability can be accepted for any damage caused by non-compliance with these instructions or any other improper use or mishandling of appliance.

## **IMPORTANT**

If you have any concerns regarding your health, consult your doctor before using this appliance. Individuals with pacemakers should not use this product. If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.



**WARNING** Only use the Battery Charger that was supplied with this appliance. Model No. GQ24-126150-AA. This Battery Charger is not designed to charge any other battery type other than Li-ion.

Never try to charge non-rechargeable batteries as they may explode or catch on fire.

**Do Not** use whilst charging.

Consumers are cautioned against trying to replace the battery. Using and/or purchasing counterfeit and/or inferior quality batteries have the potential to cause fire, personal injury, and/or damage to property. These batteries are typically found on online-auction websites or from other retailers, and are known to pose a potential safety risk.

- Batteries may explode if disposed of in fire.
- The battery is not replaceable. Do Not open the Appliance.
- Do not use or store in temperatures above 40°C



At the end of the appliance's life, it must be disposed of correctly.

Contact your local waste management authority for more information on Battery Recycling and disposal regulations in your area.

Do not dispose of electrical appliances or batteries as unsorted municipal waste, use separate collection facilities. Contact your local council for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. Please contact your local council or regional authority for more information on collection, reuse, and recycling programmes.

## Read and Save These Instructions

### 3 Introduction


Congratulations on your purchase of the Rechargeable Neck and Shoulder Massager.

It combines heat therapy together with 8 massage fingers which rotate in unison to provide you with a deep-penetrating warm and rejuvenating massage at home or away. Your body will feel more relaxed and invigorated.

The heating function can also increase the sense of well-being when used in conjunction with the rotating fingers.

## 4 Components

Heated Shiatsu heads

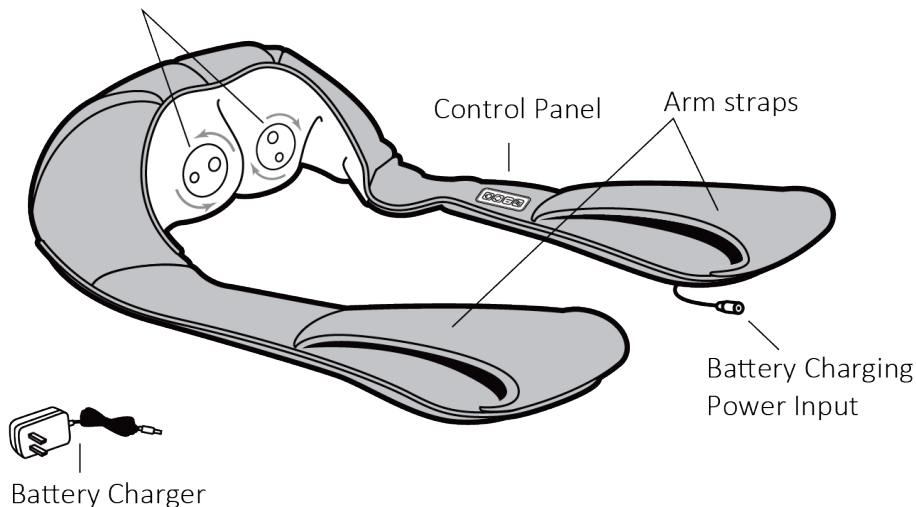

## 5 Operation Procedure

**NOTE:** If you forget to turn the massager off after the recommended maximum 15 minutes of use, the unit will automatically turn off after approximately 15 minutes.

**Do Not** use for more than 15 minutes at a time without a 30-minute cool down period in between uses otherwise, the massager will overheat.

**Do Not** use more than 2 times a day. (total of 30 minutes per day)

Always remember to turn off and unplug when not in use.

### 5.1 Charging


**DO NOT** use the massager whilst charging. Charge the battery first then use.

**Note:** Charge the Massager for at least 4 hours prior to first use.

There are no battery charge status indicator lights, so if the massager runs out of battery power, ensure to recharge for a minimum of 4 hours.

Only use the Battery Charger which was supplied with this appliance.

Take the small plug on the end of the battery charger cable and plug it into the charging port in the end of massager charging cable.

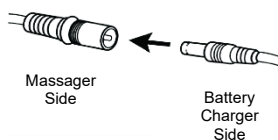

When the battery is fully charged, unplug the charger from the massager and switch off and unplug the battery charger from the power supply.


## 5.2 Massage Control Operation

### Power Button

Press to turn on the massager and activate shiatsu clockwise rotation at high level, Press again to turn off the massager

### Intensity Button

When the massager is turned on, Press to set the shiatsu intensity at high or low level

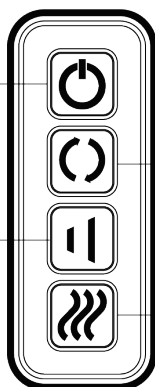

### Direction Button

When the massager is turned on, Press to shift Shiatsu clockwise rotation or Shiatsu anti-clockwise rotation

### Heat Button

When the massager is turned on, Press to turn the heat on or off


**NOTE:** Do Not use for more than 15 minutes at a time without a 30-minute cool down period in between uses.

If you forget to turn the massager off after the recommended maximum 15 minutes of use, the unit will automatically turn off after approximately 15 minutes.


**Note:** Avoid using the massager before going to sleep or in conditions where you may fall asleep when using the massager.




## 6 How to Position your Massager

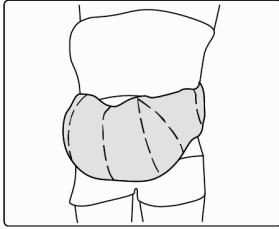


**Abdomen**

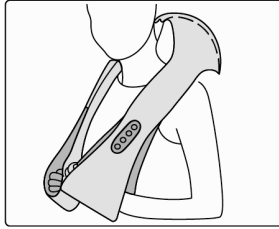


**Neck&Shoulder**

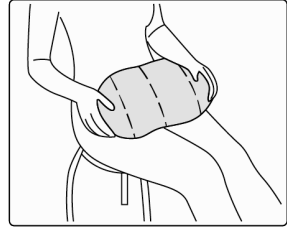


**Thigh**

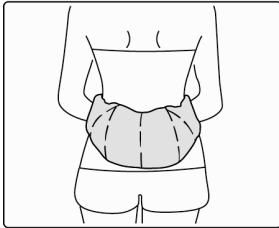

**Waist**

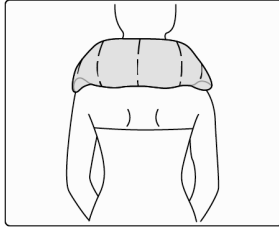

**Neck**

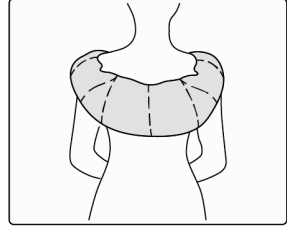

**Back**

### **IMPORTANT:**

Remove the Battery Charger from the power outlet when not in use.

For your own safety, the massager is fitted with protection against over-heating.

Heat is created in the motors when placed under great stress and excessive pressure.

Before this could create any risk, and after a runtime of 15 minutes, the massager will switch off automatically and will only allow you to switch on again after an appropriate cool down time of up to 1 hour.

## 7 Cleaning the Massager

### **CAUTION: DO NOT IMMERSE IN WATER OR LIQUID.**

- Turn off and unplug from battery charger prior to cleaning.
- Do Not immerse the Massager or battery charger in water, or any other liquid.
- Spot Clean only with a soft, very slight damp sponge.
- **NEVER** allow water or any other liquids to come into contact with the massager especially the charging socket or and battery charger.
- Never use abrasive cleaners, solvents, steel wool, scouring pads, brushes, or strong chemicals that might or might not be flammable and/or damaging to the product to clean.

**DO NOT** place in a washing machine or dishwasher.

## 8 Completion and Storage

- After you have finished using the massager, clean it as described under the Cleaning section then store safely in its box and in a dry place.
- Do Not wind the battery charger cord around the massager or charger body as this will damage the cord over time. Simply loop the cable and place next to the charger.

## 9 Technical Specifications

**Weight:** Approx. 1.8Kg Gross

**Voltage:** Input: AC 100-240V~ 50-60Hz

**Battery Charger Power:** Output: 12.6VDC 1500mA max. 18.9 Watt

**Automatic runtime:** 15 minutes

**Battery Max. runtime:** 60 minutes

**Battery Full Charging:** up to 4 hours

**Battery Pack:** 11.1Vd.c. 2200mAh Li-ion (x 3 Cells)

# 10 Troubleshooting

Problem	Cause	Remedy
No Power.	Battery has not been charged.	Plug in battery charger to massager and power point and allow to charge for a minimum of 4 hours
No Massage	Appliance not switched on	Press the On / Off button
		Refer to Section "No Power" for Remedy
Massage heads rotating at reduced speed.	The load on the massage rotating fingers is excessive.	Reduce pressure on the rotating fingers and try again.
Appliance turns off during use.	The appliance has an auto off function which activates after approx. 15 minutes of use.	Allow the product to cool for at least 30 minutes and then turn on again.
	Over-heat protection function activated. Device is switched off for your safety.	Allow the product to cool for at least 60 minutes and then turn on again.
Jerky motion of massage heads/Vibrations from device.	The load on the massage rotating fingers is excessive.	Reduce pressure on the rotating fingers and try again.
Massager turns off when Charging and in use.	You should not be using the Massager whilst it is charging.	Unplug the battery charger and allow the unit to rest for 30 minutes before recharging. <b>Do Not</b> use whilst charging.
No Indicator lights showing battery charge status	There are no indicator lights showing level of charge.	Ensure the battery is charged for at least four hours from a flat battery.

# 11 Warranty Against Defect

## 12 Month Warranty

Thank you for your purchase from Kmart.

Kmart Australia Ltd warrants your new product to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

Kmart will provide you with your choice of a refund, repair or exchange (where possible) for this product if it becomes defective within the warranty period. Kmart will bear the reasonable expense of claiming the warranty. This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or alternatively, via Customer Help at [Kmart.com.au](http://Kmart.com.au) for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to our Customer Service Centre at 690 Springvale Rd, Mulgrave Vic 3170.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

### **IMPORTANT!**

**For all technical queries or difficulties in operating the product and for spare parts, contact H.E. Group customer service 1300 105 888 (Australia) and 09 8870 447 (New Zealand).**